

Cobra SMARTsense Spirometer

12936-01

PHYWE Systeme GmbH & Co. KG Robert-Bosch-Breite 10 D–37079 Göttingen

Telefon +49 (0) 551 604-0 Fax +49 (0) 551 604-107 E-mail info@phywe.de Internet www.phywe.com



Operating instructions

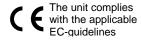


Fig. 1: 12936-01 Cobra SMARTsense Spirometer

TABLE OF CONTENTS

- 1 SAFETY PRECAUTIONS
- 2 PURPOSE AND CHARACTERISTICS
- 3 FUNCTIONAL AND OPERATING ELEMENTS
- 4 NOTES ON OPERATION
- 5 HANDLING
- **6 TECHNICAL DATA**
- 7 SCOPE OF DELIVERY
- 8 ACCESSORIES
- 9 CONFORMITY
- 10 DISPOSAL

1 SAFETY PRECAUTIONS



Caution!

- Carefully read these operating instructions completely before operating this instrument. This is necessary to avoid damage to it, as well as for user-safety.
- Only use the instrument for the purpose for which it was designed.
- Only use the instrument in dry rooms in which there is no risk of explosion.
- Protect the instrument from dust, moisture and vapours.
 Use a slightly moist lint-free cloth to clean the instrument.
 Do not use aggressive cleaning agents or solvents.
- Take care that no liquid penetrates in through the housing openings, as such penetration would result in damage to Sensor.
- Do not open the unit.

2 PURPOSE AND CHARACTERISTICS

The sensor is used for measuring the respiratory volume and volumetric flow rate and for transferring the values to a terminal device, e.g. a tablet computer, smartphone, etc., via Bluetooth or USB.

The sensor is neither suitable nor approved for medical purposes. It is intended solely for didactic purposes. This device must not be used for defined measurements on humans in order to diagnose an illness or disease, i.e. it is not intended for monitoring, treating or alleviating illnesses or diseases. As a consequence, it is not subject to the strict constraints of the Medical Devices Act.

3 FUNCTIONAL AND OPERATING ELEMENTS

3.1 Operating elements

The sensor has an on-button and two LEDs for indicating the Bluetooth and battery charge status.

On-button 🔱

To switch the sensor on and off in Bluetooth mode, the power button must be pressed for longer than 3s. If the sensor is to be connected via USB, it is not necessary to press the power button.

Bluetooth-LED *

Flashing red every 2 seconds	Not connected
Flashing green every 2 seconds	Connected to the ter- minal device
Flashing green every 4 seconds	Running measurement

Battery charge LED 🖅

Flashing red every 2 seconds	Low battery	
Illuminated red	Active charging process Charging process completed	
Illuminated red		

3.2 USB port

The battery, which is permanently installed in the sensor, is charged via the type C USB port. Furthermore, communication with a computer takes place via this interface.

4 NOTES ON OPERATION

This device fulfils all of the technical requirements that are compiled in current EC guidelines. The characteristics of this product qualify it for the CE mark.

The individual connecting leads are each not to be longer than 2 m.

The instrument can be so influenced by electrostatic charges and other electromagnetic phenomena (HF, bursts, indirect lightning discharges) that it no longer works within the given specifications. Carry out the following measures to reduce or eliminate the effect of such disturbance: Ensure potential equalization at the PC (especially with Laptops). Use screening. Do not operate high frequency emitters (e.g. radio equipment or mobile radiotelephones) in the immediate vicinity. When a total failure of the instrument occurs, unplug it and plug it back in again for a reset.

5 HANDLING

This section describes the start-up of the sensor and the recording of measurement data. Please read this section thoroughly in order to avoid failures or operating errors.

5.1 Charging process

Use a USB-C cable to connect the sensor to a computer or USB charger (not included).

During the charging process, the battery charge LED lights up red. When the charging process is complete, the battery charge LED lights up green. The charging time for a completely discharged battery is 3 hours maximum.



Disconnect the charger at the latest four hours after the completion of the charging process. Otherwise, the service life of the battery may be negatively affected.

5.2 Start-up

Attach the supplied mouthpiece with filter to the sensor.



Fig. 2: Connecting the mouthpiece with filter

Switch on the sensor by pressing the power button for more than 3s. Now the Bluetooth LED flashes red. Start the software and select the sensor.

If the sensor is to be used via the USB interface, it does not need to be switched on. The sensor is connected directly to the end device using the supplied USB cable.

There is a 9-digit code on the back of the sensor (Fig.3). The last 4 digits of the code are displayed as the sensor name in the software (Fig.4). This enables the precise assignment of the sensors within the software.



Fig. 3



Fig. 4

Selection of the sensor via the Bluetooth interface

Make sure that the Bluetooth interface is activated on the terminal device (PC/Tablet/Smartphone) and that the software is allowed to access the interface.

After the sensor has been selected in the software, the LED flashes green to indicate that the connection has been established correctly. After the sensor has been coupled with the software, the sensor is no longer visible to other users in the software, and therefore can no longer be selected.

If the sensor is switched on and not connected, it switches off automatically after 5 minutes.

Selection of the sensor via the USB interface

For this purpose the sensor must be plugged into the USB port of the end device. It is not necessary to switch on the sensor. The sensor is automatically recognized and displayed. It can be selected and connected directly.

5.3 Recording of measurement data

Measuring principle:

Inside the tube, the air flow hits several small lamella. The result is a slight flow resistance, leading to a small difference in pressure. This pressure difference is evaluated electronically. In accordance with the Hagen–Poiseuille law, the pressure drop upstream and downstream of the lamella in laminar flow is directly proportional to the volume flow Q (flow rate).

The volume V is calculated by integrating the flow rate as a function of time.

Set flow rate to zero:

Depending on how you hold the sensor, an offset value of the flow rate may be displayed. This can cause the calculated volume to increase continuously even though you are not breathing into the sensor. To prevent this, set the flow rate to zero using the software.

Measurement:

Stand upright and start the measurement. Exhale and inhale evenly through the attached cardboard tube. Then, inhale as deeply as possible and then exhale completely, if possible. Stop the measurement. This spirogram provides information about the various parameters of the respiratory volume, e.g. the vital capacity and tidal volume.

6 TECHNICAL DATA

Operating temperature range: 5 - 40°C Rel. humidity < 80%

Flow Rate

Measuring range	-10 l/s ~ 10 l/s
Resolution	0.01 l/s
Accuracy	±3%

Volume:

Measuring range	-15 ~ 15
Resolution	0.01
Accuracy	±3%
	1
Max. date rate	200 Hz
Battery capacitiy	250 mAh
Max. wireless range (open field)	30 m
Dimensions (length x width x height)	120x 40 x 37 mm
Weight	63 g

7 SCOPE OF DELIVERY

The extent of delivery is as follows

•	Cobra SMARTsense Spirometer	12936-01
•	USB cable type C	07935-00
•	1x Mouthpiece including filter	12936-20

Operating instructions

8 ACCESSORIES

The following accessories are available:

•	Cobra SMARTlink	12999-99
•	USB-charger	07934-99
•	USB cable type C	07935-00
•	USB-Bluetooth-Adapter	07936-00
•	Software measureLAB	14580-61
•	Mouthpiece including filter	12936-20
•	Cardboard mouthpiece with filter paper	12936-10

• Free measureApp available from supplier portals



9 CONFORMITY



PHYWE Systeme GmbH & Co.KG hereby declares that the radio system type 12936-01 complies with the 2014/53/EU directive. The complete text of the EC Declaration of Conformity is available at the following Internet address:

www.phywe.com/en/ec-declaration

10 DISPOSAL

The packaging mainly consists of environmentally-friendly materials that should be returned to the local recycling stations.



Do not dispose of this product with normal household waste. If this unit needs to be disposed of, please return it to the address that is stated below for proper disposal

PHYWE Systeme GmbH & Co. KG Abteilung Kundendienst Robert-Bosch-Breite 10 D–37079 Göttingen

Telefon +49 (0) 551 604-274 Fax +49 (0) 551 604-246